

DiamondCare® Outlet Surface Mount Enclosures

SPECIFICATION

Surface Mount Enclosures

Surface mount enclosures permit installation of DiamondCare® medical gas outlets and slides on a wall surface where a recessed wall outlet cannot be installed.

Surface Mount Gas Outlet Enclosures

The Surface Mount Gas Outlet Enclosure is designed to mate with a DiamondCare® Console Outlet for any available gas. The DiamondCare® Console Outlet finish and rough-in (standard piped) assemblies must be ordered separately.

The Surface Mount Gas Outlet Enclosure consists of a powder-coated, light-neutral colored, metal enclosure, a mounting bracket and two screws. The front of the enclosure kit includes a cut-out for mounting the desired DiamondCare® Console frontbody. The mounting bracket is used to affix the DiamondCare® Console Outlet rough-in assembly (standard piped) to the wall. The enclosure is held in place by attaching the console frontbody to the rough-in assembly.

Surface Mount Slide Enclosures

The Surface Mount Slide Enclosure provides a slide attachment that consists of a powder-coated light neutral, metal enclosure, a mounting bracket and slide with two mounting screws.


Features


- Permits wall mounting of DiamondCare® medical gas services and slides where a recessed wall outlet cannot be installed.
- Can be used with any available DiamondCare® gas service and key style outlet.

Enclosure	Product Number
DiamondCare® Wall Outlet - Surface Mount Gas Outlet Enclosure	❑ 4107 2011 80
DiamondCare® Wall Outlet - Surface Mount Slide Enclosure	❑ 4107 2011 81

NOTE:

A DiamondCare® Console Outlet must be ordered separately for use with the Surface Mount Gas Outlet Enclosure.

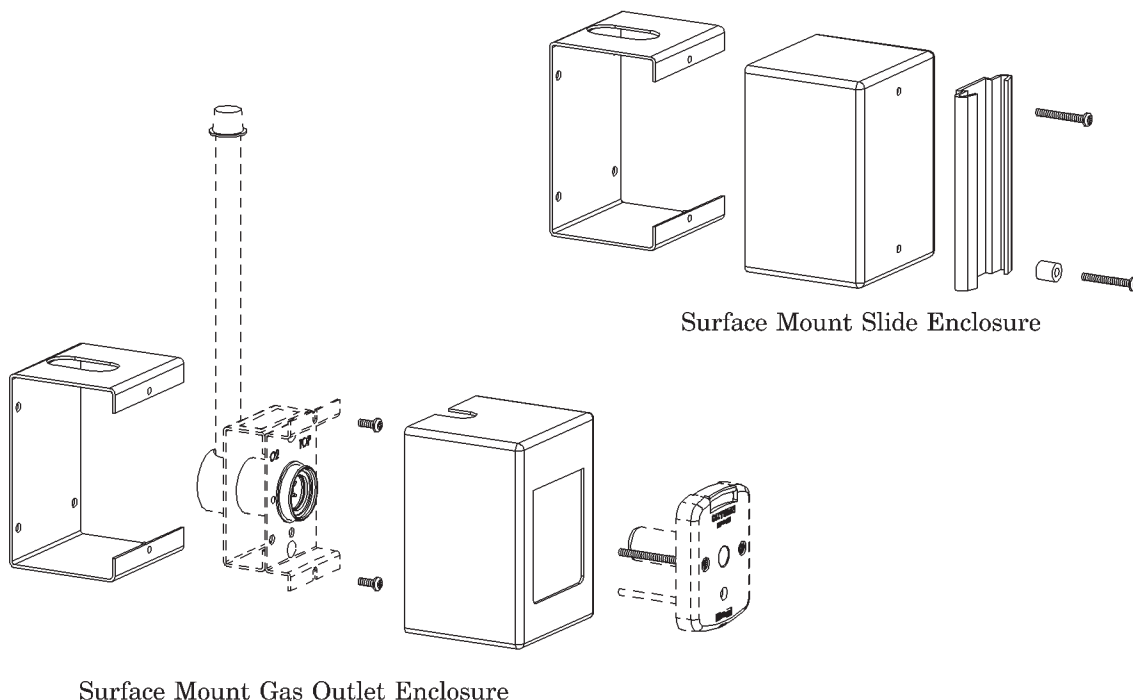

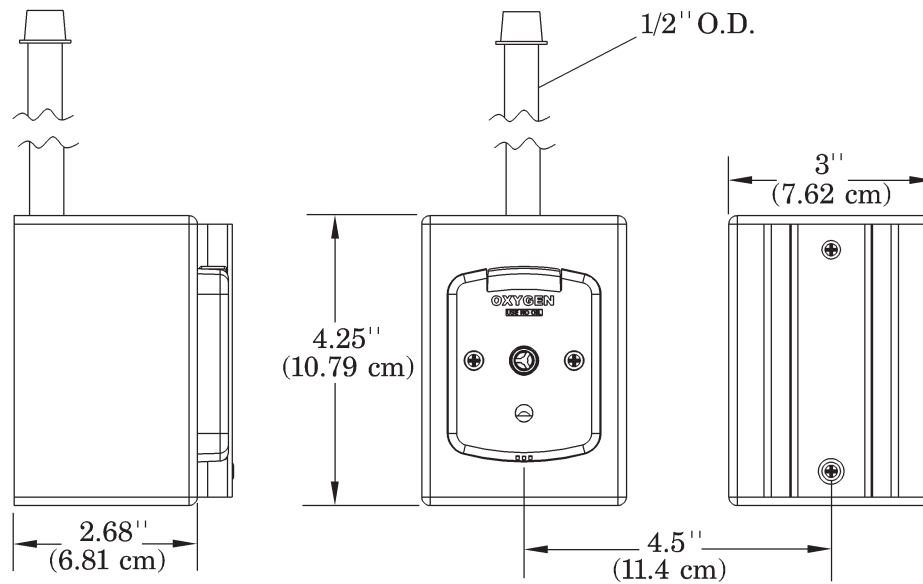Product Dimensions

For more information, contact BeaconMedaes Customer Service at 1-888-4MEDGAS (463-3427).

t011a002